


# AddOn®

hydrophilic • refractive

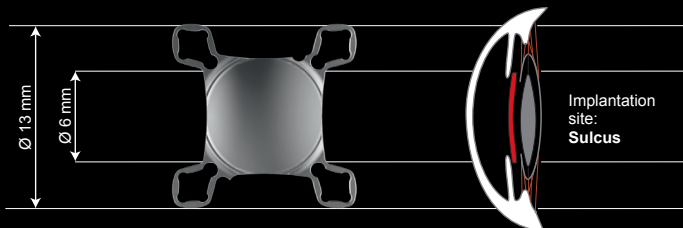

## Q functional

ONV function\*

IOL function

Distant function

Basis function

## Features

One piece, foldable intraocular lens (IOL), biocompatible hydrophilic acrylic (25% water content), as additional IOL in combination with a Basis IOL.

- |                      |                           |
|----------------------|---------------------------|
| Convex/Concave Optic | ▶ Cell Stop/Non IOL Touch |
| 4 Flex-Haptics       | ▶ Sulcus Fit              |
| Non Torque Design    | ▶ Rotational Stability    |
| Square Design        | ▶ Non Iris Capture        |

Optical power:



AddOn®	Sph.	Steps
A4SW00	-10.0 D to +10.0 D	0.25 D

\*ONV functions: The five functions of the human eye which are required to achieve optimal natural vision (ONV) (see rear page)